Addendum

223 Sample Changer User's Guides

Part Numbers LT1931 and LT1932

This document contains corrections and additions to information provided in the *223 Sample Changer User's Guides*.

- **Safety** and **Sécurité** sections are added.
- **Fuse Replacement** replaces the **Replacing a Fuse** section in the **Maintenance** chapter.
- **Fuse Installation** replaces the **Fuses** section in the **Installation** chapter.
- **Fuse Part Numbers** replaces the fuse specifications in the **Replacement Parts** section of the **Replacement Parts and Accessories** chapter.

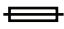
Safety

The instrument is intended to be used in a laboratory by trained technical personnel. For safe and correct use, both operating and service personnel should follow the instructions contained in the user's guide when installing, cleaning, and maintaining the instrument.

The following safety precautions must be observed during all phases of operation, service, and repair of the instrument. Failure to comply with these precautions or with specific warnings elsewhere in the user's guide violates safety standards of design, manufacture, and intended use of the instrument. Gilson assumes no liability for the customer's failure to comply with these requirements.

The instrument has been certified to safety standards required in Canada, Europe, and the United States. Refer to the rear panel label on the instrument or the Declaration of Conformity document for the current standards to which the instrument has been found compliant.

The following electronic and hazard symbols may appear on or inside the instrument.

Symbol	Explanation	Symbol	Explanation
~	Alternating current Courant alternatif Wechselstrom		Direct current Courant continu Gleichstrom
	Protective conductor terminal Borne de terre de protection Schutzleiteranschluss		Caution Attention Vorsicht
	Electrical power ON Sous tension Netzschalter ein		Caution, risk of electric shock Attention, risque de choc électrique Vorsicht, Elektroschockgefahr
O	Electrical power OFF Hors tension Netzschalter aus		Fuse Fusible Sicherung
 KEEP HANDS CLEAR OF PROBE!	Keep hands clear of probe Garder les mains éloignées de l'aiguille Halten Sie Hände fern von der Nadel		

Voltage

Access to the rear panel is necessary. The instrument must be detached from all voltage sources before service, repair, or exchange of parts. For normal operation, the instrument is to be grounded through the AC line cord and power supply provided. Failure to do so can result in a potential shock hazard that could result in serious personal injury.

Use only fuses with the rated current and of the specified type as listed on the rear panel label on the instrument. The instrument must only be operated with the voltage specified on the rear panel label of the instrument and with the grounded AC line cord and power supply provided.

Solvents

Observe safe laboratory practices when handling solvents. If working with hazardous solvents or flammable liquids, ensure that there is proper ventilation and that adequate protection such as safety glasses, gloves, and protective clothing are used.

If dangerous liquids are used, precautions should be taken to limit potential hazards from leaks and/or spillage through the use of a non-flammable tray or use of a fume hood, etc.

If there is the potential of explosive gases being developed, a fume hood or other means should be used to safely manage that risk.

Refer to the Material Safety Data Sheets for the solvents before use.

Replacement Parts

Be sure to use only replacement parts mentioned in the user's guide. Do not repair the instrument or change parts not listed in the user's guide. If it is necessary to replace parts not listed, please contact your local Gilson representative.

Probes

Keep hands clear of probe to avoid risk of personal injury by piercing while operating the instrument. Because the probe installed on the Z-arm may contain dangerous substances, do not interfere in the work area of the instrument until it has completed its procedures.

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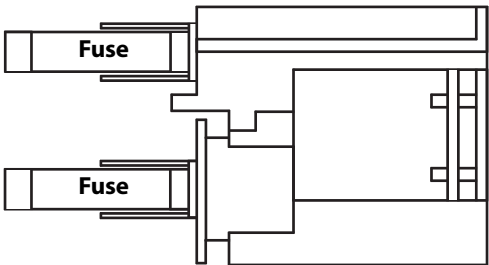
Fuse Installation

- 1 Locate the fuse drawer and two of the supplied fuses.
- 2 Insert the fuses in the fuse drawer.
- 3 Insert the fuse drawer into its receptacle on the rear panel.

Fuse Replacement

- 1 Power off the instrument and disconnect the power cord.
- 2 Locate or order replacement fuses. (Extras were provided with the instrument.)
- 3 Place a small screwdriver or a fingernail under the tab on the fuse drawer to detach it.
- 4 Remove the fuse drawer from its receptacle on the rear panel.
- 5 Replace both fuses. Use only fuses with the rated current and specified type as listed below and on the rear panel of the instrument.
- 6 Insert the fuse drawer into its receptacle on the rear panel.

Fuse Part Numbers

<i>Part Number</i>	<i>Description</i>
6730204007	FUSE, 2.0 AMP T-2 SLO-BLO
6770100411	FUSE DRAWER